





Biotek is a global medical device company with a strong focus on new product development in the most advanced orthopedic sector: SPORTS MEDICINE

Biotek continues to experience unprecedented growth and demand for our products throughout the world; we currently have presence in all six continents around the world. However, Biotek remains a privately held company with a family business culture committed to delivering uncompromising quality to the health care professionals who use our products, and ultimately, the millions of patients whose lives we impact.

year founded
2003



40+
countries
our products are
used in




100+
Sales partners
globally



MINI-VIM®

Ligament Anchor
with Needle (Titanium)

ELBOW
SPORTS MEDICINE

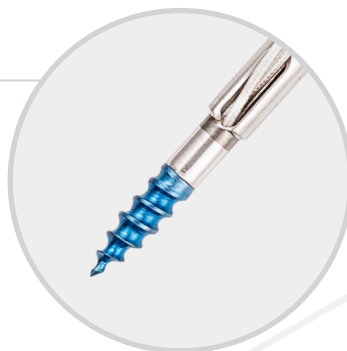

MINI-VIM® Ligament Anchor with Needle (Titanium)

- The MINI-VIM® Ligament anchors are titanium anchors available in 2.8mm, 3.5mm diameter and come pre-loaded with USP #2 suture and needle.
- The implant can be manually inserted using the handled version. A drill is available for use in hard bone if desired.
- **Anchor Dimensions:** 2.8mm × 11.5mm, 3.5mm × 11.5mm
Type: Screw in suture anchor
Suture Colour: White

ORDERING INFORMATION

BAS-9040.28F	MINI-VIM® Ligament Anchor, Titanium loaded with 1 pc Fiber with needle, Dia. 2.8mm
BAS-9083.35F	MINI-VIM® Ligament Anchor, Titanium, loaded with 1 pc Fiber with needle, Dia. 3.5mm

RELATED INSTRUMENTS

BAS-9026	Drill bit, Dia. 2.2mm for MINI-VIM® Ligament Anchor BAS-9001.28F & BAS-9040.28F
BAS-9114	Drill bit, Dia. 2.8mm for MINI-VIM® Ligament Anchor BAS-9031.35F & BAS-9083.35F

MINI-VIM PK®

Ligament Anchor
(PEEK, Knotless)

ELBOW
SPORTS MEDICINE



MINI-VIM PK® Ligament Anchor (PEEK, Knotless)

- The MINI-VIM PK® Ligament anchor, Knotless are available in diameter of 2.8mm & 3.5mm.
- It is a non-absorbable PEEK suture anchor with a material eyelet which provides superior abrasion resistance due to PEEK's low coefficient friction.
- The knotless technology of MINI-VIM PK® Ligament anchor provides a step-saving alternative to conventional "knotted" suture anchors and eliminates "knot stacks" associated with soft tissue irritation.
- **Anchor Dimensions:** 2.8mm x 12.0mm, 3.5mm x 12.0mm
Type: Tap-in/Bang-on suture anchor
Suture Colour: White-Green (Holding Suture)

ORDERING INFORMATION

BAS-9087.28F	MINI-VIM PK® Ligament Anchor, PEEK, Knotless, Dia. 2.8mm
BAS-9192.35F	MINI-VIM PK® Ligament Anchor, PEEK, Knotless, Dia. 3.5mm


RELATED INSTRUMENTS

BAS-9111	Drill bit, Dia.2.5mm for MINI-VIM PK® Ligament Anchor, Knotless-BAS-9087.28F
BAS-9127	Sleeve for MINI-VIM PK® Ligament Anchor, Knotless-BAS-9087.28F

FIBERKNOT®

Ligament Anchor
with Needle (Soft)

ELBOW
SPORTS MEDICINE

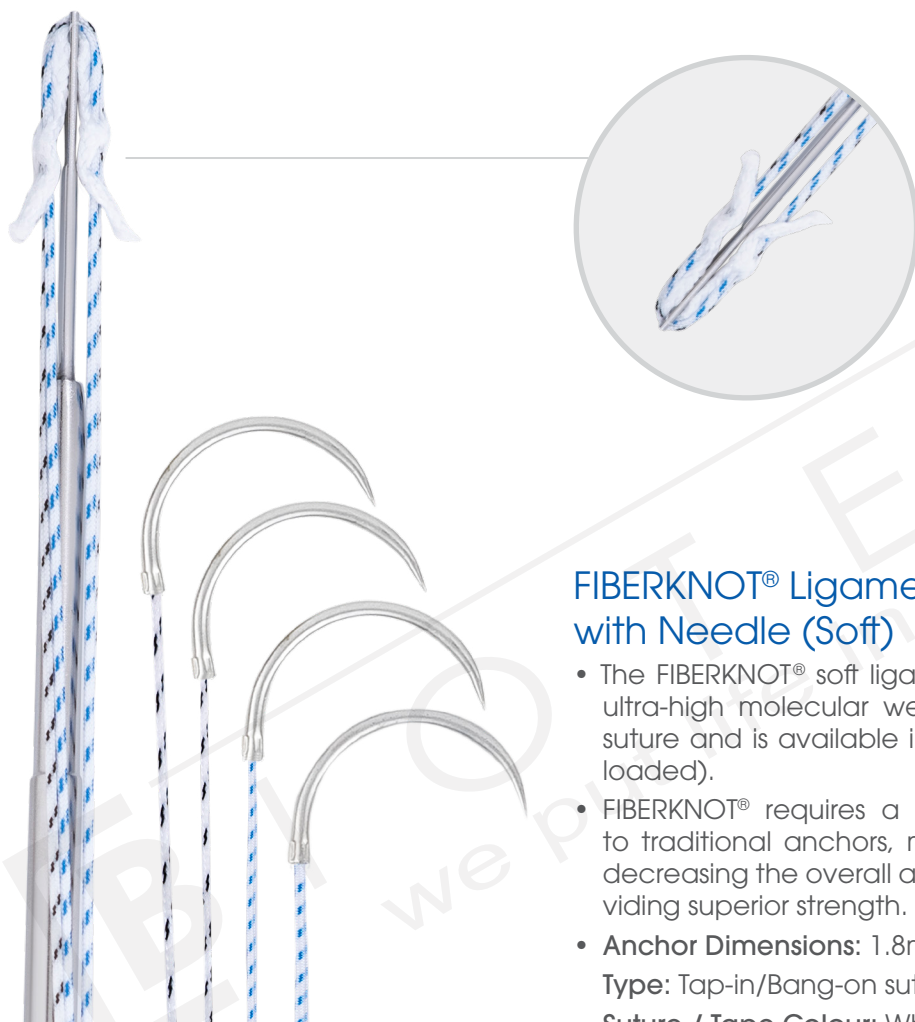

FIBERKNOT® Ligament Anchors with Needle (Soft)

- The FIBERKNOT® soft ligament anchor is composed of ultra-high molecular weight polyethylene (UHMWPE) suture and is available in diameter of 1.8mm (double loaded).
- FIBERKNOT® requires a smaller pilot hole compared to traditional anchors, minimizing bone removal and decreasing the overall anchor footprint, while still providing superior strength.
- **Anchor Dimensions:** 1.8mm × 15.0mm
Type: Tap-in/Bang-on suture anchor
Suture / Tape Colour: White-Blue, White-Black

ORDERING INFORMATION

- | | |
|-------|--|
| 10095 | FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber with needle Dia. 1.8mm |
| 10030 | FIBERKNOT® Soft Ligament Anchor loaded with 1 pc White-Black Fiber and 1 pc White-Blue Fiber Tape size 1.4, With Needle Dia. 1.8mm |

RELATED INSTRUMENTS

- | | |
|-------|---|
| 10066 | Drill Bit Dia 1.6mm , Len.- 100 mm for FIBERKNOT® Soft Ligament Anchor Dia 1.5mm and 1.8mm with Needles |
|-------|---|



AK-FIX Fixation Button


- The AK-Fix Fixation Button is a low profile Titanium button with slots enabling the insertion of the Biofiber loop.
- The oblong shaped button can be used with USP size 5 BioFiber Loop to treat Distal Biceps Repair.
- Button Dimensions: 10.0mm x 3.0mm x 1.3mm

ORDERING INFORMATION

10104	AK-FIX Fixation Button
10122	AK-FIX Fixation Button, with Inserter

RELATED INSTRUMENTS

10144	Suture Passing Inserter for AK-FIX Fixation Button
-------	--

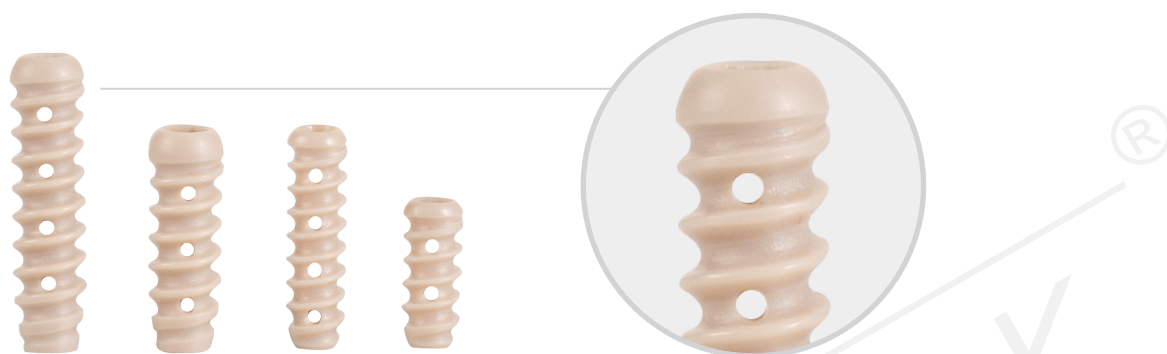

SOFTFIX-PK® Interference Screws, PEEK

- Non-absorbable SOFTFIX-PK® interference screw incorporates PEEK-VESTAKEEP® polymer from Evonik® which exhibits an ideal combination of strength, stiffness and toughness that surpasses absorbable materials.
- It provides the interference screw fixation, while eliminating the need of trans osseous tunnels in tendon repairs and ligament reconstruction.
- It can be used for Distal Biceps Repair.
- Dimensions: Dia. 4.0mm x Len. 10.0mm & 15.0mm and Dia. 5.0mm x Len. 15.0mm & 20.0mm

ORDERING INFORMATION

10053	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 4.0mm, Len. 10mm
10054	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 4.0mm, Len. 15mm
BAK-7164.15	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 5.0mm, Len. 15mm
BAK-7164.20	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 5.0mm, Len. 20mm

RELATED INSTRUMENTS

10008	Screwdriver for SOFTFIX-PK® Interference Screw, Cannulated, Dia. 4.0mm
BAK-7223	Screwdriver for SOFTFIX-PK®, Dia. 5.0mm Interference Screws, Cannulated
10085	Fluted Reamer Cannulated, Dia. 4.0mm
10086	Fluted Reamer Cannulated, Dia. 5.0mm
BAL-024	Threaded Guidewire Dia. 1.2; Len. 180mm, SS
10094	Guide Wire, Blunt Tip for Dia. 4.0mm SOFTFIX-PK® Interference Screws Dia. 0.8mm; Len. 200mm (SS)
10096	Guide Wire, Blunt Tip for Dia. 5.0mm-8.0mm SOFTFIX-PK® Interference Screws Dia. 1.2mm; Len. 200mm (NITINOL)
10065	Drill Tip Suture Passing Pin, Dia. 2.4mm; Len. 200 mm



BioFiber®

BioFiber® is a braided sterile surgical fiber prepared from fibres of ultra-high molecular weight long chain linear polyethylene. BioFiber® elicit a minimal of inflammatory reaction in tissues, followed by gradual encapsulation of the fiber by the fibrous tissue.

ORDERING INFORMATION

BF-2W	BioFiber® Fiber USP Size 2 with MO-6 needle White, Len. 36"	BF-WWB	BioFiber® Fiber USP Size 2 Len. 36" Pack contains 2 units (White with needle & White-Blue with needle)
BF-2WC	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Blue, Len. 36"	10068	BioFiber® Fiber USP Size 2-0 White, Len. 75cm
BF-2WK	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Black, Len. 36"	10069	BioFiber® Fiber USP Size 0 White, Len. 75cm
BF-2WCK	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Blue-Blue-Black, Len. 36"		

BioFiber® Tape



BioFiber® Tape

BioFiber® Tape is an ultra-high strength 2 mm width tape using a similar long chain polyethylene structure as the BioFiber® suture. The broad footprint of the BioFiber® Tape is appropriate for repairs in degenerative cuff tissue where tissue pull-through may be a concern.

ORDERING INFORMATION

BF-2BST	BioFiber® Fiber Tape, size 2 White-Black, Len. 39"
BF-2ST	BioFiber® Fiber Tape, size 2, White-Blue, Len. 39"

BioFiber® Loop



BioFiber® Loop

BioFiber® Loop is a continuous loop of BioFiber® with straight needle. The needle is easy to handle and move freely on the suture to recenter itself after passing through tissue and facilitating even tension. Using BioFiber® Loop in conjunction with the Whipstitch technique drastically reduces the time spent and improves the fixation strength.

ORDERING INFORMATION

BF-2WSB	BioFiber® Fiber Loop, Size 2 with straight needle, Len. 6cm White-Blue, Len. 24"
---------	--






BIOTEK - CHETAN MEDITECH PVT. LTD.

MD-04, Sanand GIDC-II,
Ahmedabad, Gujarat-382110. India.

Email: contact@biotekortho.com

Follow us on:   



DISTRIBUTED BY