


# Biceps Tenodesis Repair

## VIMFIX<sup>®</sup>-BT Ligament Anchor (Knotless)

The VIMFIX<sup>®</sup>-BT Tenodesis implants are designed for all-arthroscopic proximal biceps tenodesis.

The efficient system is designed to save steps and minimize the length of the procedure. The implants feature a unique PEEK tip that is used to steer the biceps tendon to the bottom of the bone socket without the need to externalize or whipstitch the tendon. There is no need to predetermine tendon length, as position and tension are easily set and visualized prior to final anchor insertion. Strong fixation is obtained by simply advancing the preloaded VIMFIX<sup>®</sup>-BT Tenodesis Screw. The system also features Tenodesis Fluted Solid Reamers with built-in guide tips that eliminate the need to ream over a guide pin.

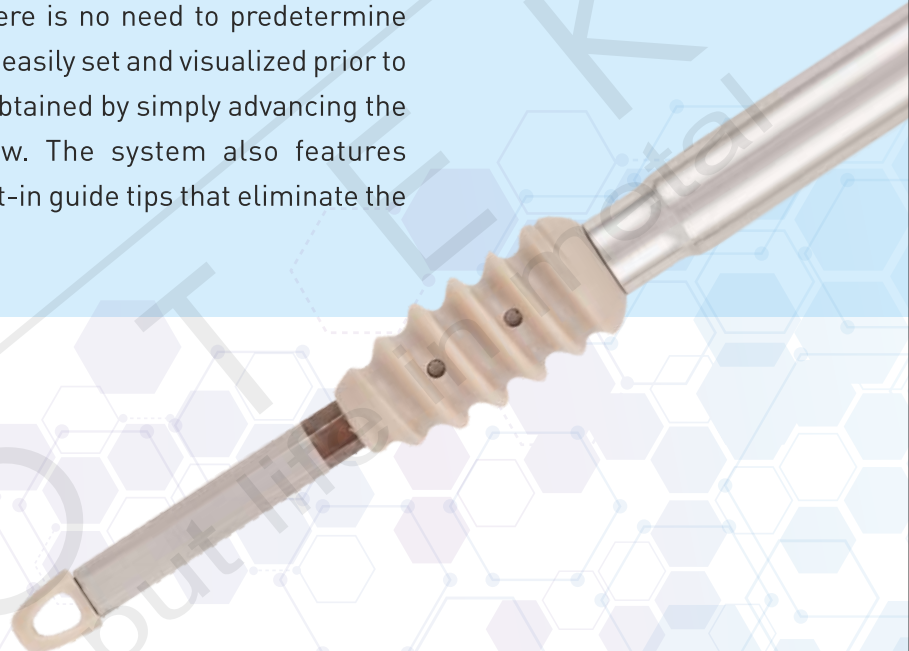

### VIMFIX<sup>®</sup>-BT Ligament Anchor

BAS-9093.50F	VIMFIX <sup>®</sup> -BT Ligament Anchor, PEEK, knotless, Dia. 5.0mm
BAS-9094.60F	VIMFIX <sup>®</sup> -BT Ligament Anchor, PEEK, knotless, Dia. 6.0mm
BAS-9095.70F	VIMFIX <sup>®</sup> -BT Ligament Anchor, PEEK, knotless, Dia. 7.0mm
BAS-9096.80F	VIMFIX <sup>®</sup> -BT Ligament Anchor, PEEK, knotless, Dia. 8.0mm

### Related Instrumentation

BAS-9113	Graft Sizer for Biceps Tendon, Dia. 5- 6-7-8 mm
BAS-9104	Fluted Reamer Solid, Dia. 5.0mm
BAS-9105	Fluted Reamer Solid, Dia. 6.0mm
BAS-9106	Fluted Reamer Solid, Dia. 7.0mm
BAS-9107	Fluted Reamer Solid, Dia. 8.0mm

